

HASO FIRST

Wavefront sensor
The Chameleon

On demand wavelength High accuracy Best cost performance ratio







HASO FIRST +

The HASO Shack-Hartmann Wavefront Sensor optimized for one wavelength, the one you really need.

The HASO FIRST is now faster and has an improved spatial resolution while keeping the same accuracy and optimized price point.



Compatible with the Optical Engineer Companion modular system: easily combine the accessories you need.

APPLICATIONS

Successfully used in the most demanding applications in optical metrology, microscopy, and laser diagnostics, the HASO FIRST performs multiple functions:

- + Quantify the aberrations of an optical system
- + Align optical systems to ensure that it works optimally
- + Predict the performance of optical systems in terms of focusing capability or imaging quality
- + Quantify the effects of temperature and gravity on system performance
- + Verify that the optics comply with specifications
- + Drive a wavefront corrector to correct for system aberrations
- + Check whether the optical mount overly distorts the optics

FEATURES

- + Beam collimation with an accuracy better than 200 m radius of curvature
- + A 20 mm focal length measurement with a sensitivity of 1 μ m RMS
- + Direct wavefront acquisition of converging and diverging F/5 beams with an accuracy of $\lambda/100$ RMS including astigmatism and high-order aberrations
- + Control and adjustment of axial laser beam deviation better than 5 $\mu rad \; RMS$
- \pm 50 nm calibration bandwidth or extended wavelength range optional: \pm 150nm around the calibration wavelength



SPECIFICATIONS*

OPERATING SPECS

Aperture dimension Number of microlenses Maximum acquisition frequency One wavelength ± 50 nm in the range Minimum power External trigger TTL signal Operating system

OPTICAL SPECS

Repeatability
Absolute wavefront measurement accuracy
Spatial sampling
Local radius of curvature dynamic range


MISC

Dimensions (Height x Width x Length) Weight for USB version Working temperature Interface Power consumption 44 x 36 125 Hz (USB 3.0) or 30 Hz (with GigE converter) 350 - 1100 nm 0.15 nW TTL signal Windows 10 & 11

< λ /200 RMS ~ λ /100 RMS ~ 100 μ m \pm 0.008 m to \pm ∞

4.5 x 3.7 mm²

42 x 47 x 60 mm³ (USB 3.0) 200 g 15 - 30 °C USB 3.0 or optional GigE converter

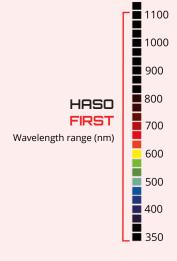



HASO FIRSTDynamic range at $\lambda = 635$ nm



 $\hbox{*Subject to changes without further notice}\\$

DIMENSIONS** (mm)

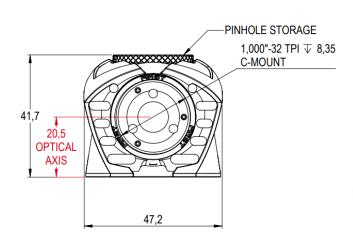

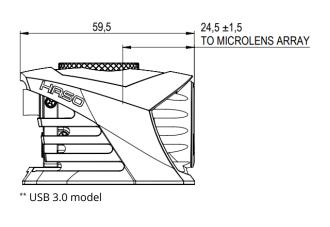

SOFTWARE

WAVEVIEW™ Metrology Software

WAVEVIEW[™] is the most advanced wavefront measurement and analysis software.

It offers more than 150 features and tools optimized for a wide range of highly demanding applications.

Options:


- + Extensions for PSF, MTF, M² and Strehl ratio
- + Optional SDK in C/C++, LabVIEW and Python


WAVETUNE™ Adaptive Optics Software

WAVETUNE™ is a unique software that seamlessly combines wavefront measurement and correction features with extensive instrument diagnostics. It is perfectly adapted to our HASO wavefront sensors, ILAO STAR, MIRAO and mu-DM deformable mirrors, as well as to a wide range of active components.

Options:

+ Optional SDK in C/C++, LabVIEW and Python

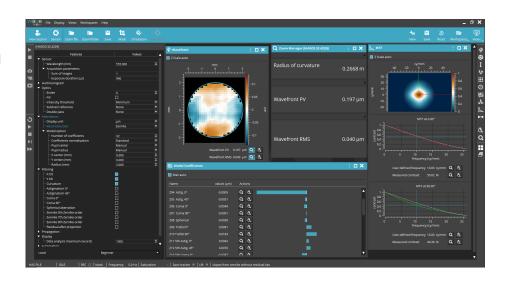

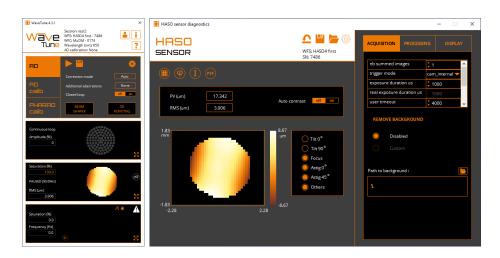

CONTACT US

Imagine Optic Headquarters
18, rue Charles de Gaulle
91400 ORSAY · France
Phone +33 (0)1 64 86 15 60
sales@imagine-optic.com
www.imagine-optic.com

